


# Chamber Safety -Temperature

Thermocouples placed at **A** & **B**

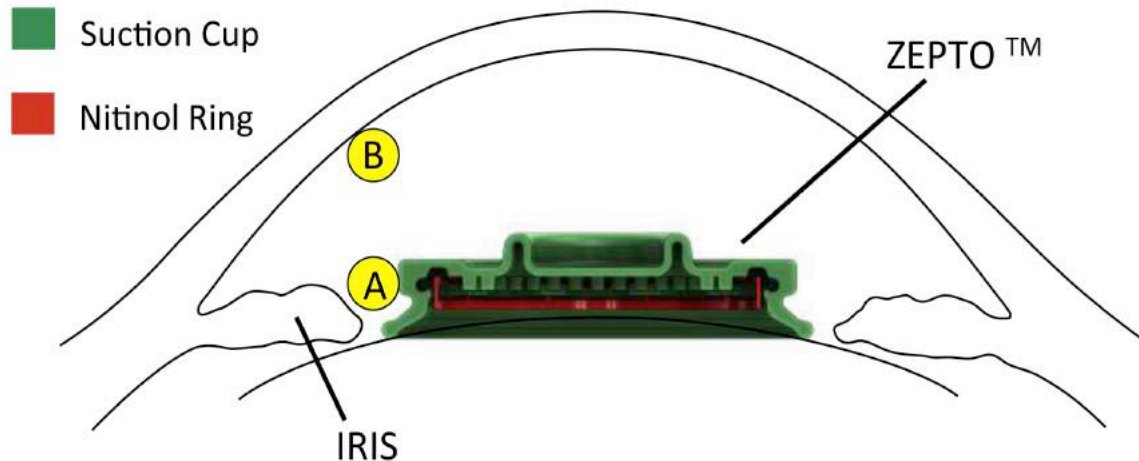

	duration; peak
<b>A</b>	2.94 secs; 1.76 <sup>0</sup> C
	3.59 secs; 1.91 <sup>0</sup> C
	3.69 secs; 2.01 <sup>0</sup> C
<b>B</b>	2.69 secs; 0.89 <sup>0</sup> C
	4.47 secs; 0.53 <sup>0</sup> C
	4.60 secs; 1.32 <sup>0</sup> C

<b>A</b>	2.94 secs; 1.76 <sup>0</sup> C
	3.59 secs; 1.91 <sup>0</sup> C
	3.69 secs; 2.01 <sup>0</sup> C
<b>B</b>	2.69 secs; 0.89 <sup>0</sup> C
	4.47 secs; 0.53 <sup>0</sup> C
	4.60 secs; 1.32 <sup>0</sup> C

No significant temperature change in eye

- Low energy capsulotomy
- Additional shielding by silicone cup